Waterpik™ Water Flosser Model WF-26 | 27 | 28 Series

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www.waterpik.com

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place or drop the product into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not handle the user supplied USB power supply (power adapter) with wet hands.

- Unplug the user supplied USB power supply (power adapter) when not in use.
- Check the charging cord and user supplied USB power supply (power adapter) for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only for its intended use as described in these instructions or as recommended by your dental professional.

- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. at 1-800-525-2774 (US) 1-888-226-3042 (Canada) or visit www.waterpik.com (US) www.waterpik.ca (Canada).
- A user supplied USB-A wall charger (power adapter), safety certified (IPX4, Class II) with a rated output of 5.0 V $\overline{=}$, 3.0 Amps Max is required to charge this product.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Keep the charging cord away from heated surfaces.
- Never drop or insert any object into any hose or opening.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Fill reservoir with water or other dental professional recommended solutions only.
- Only use tips and accessories recommended by Water Pik, Inc. as identified in this instruction manual.
- If the charging cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Do not direct water from the water flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.

- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- Do not use the water flosser for more than five minutes in each two-hour period.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
- There are no user serviceable electrical parts in this product and it does not require maintenance.
- Remove the charging cord from the unit and replace charging port cover prior to filling the reservoir.
- Ensure contacts on the cord, the unit and the power adapter are clean and dry while the unit is being charged.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet. Contact your doctor with any questions.
- Do not expose the battery or appliance to excessive temperatures.

SAVE THESE INSTRUCTIONS

INTENDED USE

Waterpik™ water flossers with PRECISIONPULSE TECHNOLOGY are proven to deliver an effective combination of water pressure and pulsations to remove plaque bacteria that can cause gingivitis, cavities, and bad breath.

Product Description	7
Getting Started	8
Cleaning and Troubleshooting Your Waterpik™ Water Flosser	11
Battery Care, Battery Disposal/Recycling, Battery Removal	12
Limited Two Year Warranty	13


Symbols Glossary (Not all symbols apply to all products)

Explanation of symbols used on Water Pik, Inc. product labels and packaging.


Applicable symbols for a specific product are found on that product's labels and packaging.

Symbol	Definition	Symbol	Definition	Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Detachable Power Supply		Do not use this device in a shower, bathtub and do not immerse in any liquid.		Separate collection for waste electrical and electronic equipment.
	Manufacturer		Direct Current (DC) Voltage		Product contains a lithium ion battery		Contains a recyclable Li-ion battery
IPX7 or IPX4	Ingress Protection Code		Alternating Current (AC) Voltage		CE marking of conformity		Contains a recyclable Ni-MH battery
	IEC Class II Equipment		Do not immerse this device in any liquid.		Regulatory Compliance Mark		

1. EasySpin™ 360 degree tip rotator
2. Power button
3. Pressure indicators
4. Power handle
5. Charge port/Charge port cover
6. Reservoir
7. Tip
8. Reservoir fill door
9. Battery charge indicator
10. Tip eject button

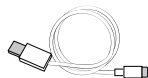Charging Cable
20033894-2 White

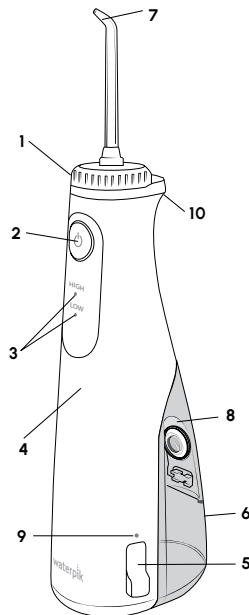*USB power adapter (not included)

TIPS

A. Orthodontic Tip* (OD-100E)

- Braces/General use

B. Pik Pocket™ Tip* (PP-100E)

- Periodontal pockets/Furcations

C. Plaque Seeker™ Tip* (PS-100E)

- Implants/Crowns/Bridges/
Retainers/General use

D. Handheld Cordless

Precision Tip (JT-450E)**

- General use






E. Tongue Cleaner* (TC-100E)

- Fresher breath


G. Travel Plug

- Keeps excess water from
draining out while traveling

NOT ALL TIPS INCLUDED IN ALL MODELS

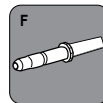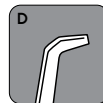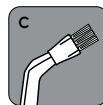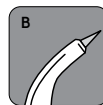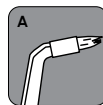*Replace Every 3 months **Replace Every 6 months
Replacement tips/attachments may be purchased online at
www.waterpik.com (US) or www.waterpik.ca (Canada), or by phone from
Water Pik, Inc. at 1-800-525-2774 (US) or 1-888-226-3042 (Canada).

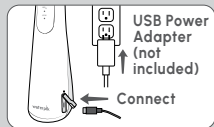Charging Unit Prior to First Use

Connect the charging cable to a USB power adapter (not included), then plug into a wall outlet. Attach the other end of the cable to the unit.

Charge for **4 hours** before first use. A full charge lasts up to 9 weeks with typical use.

Replace the charge port cover after charging and ensure it's flush with the handle.

Note: Refer to the IMPORTANT SAFEGUARDS section for required USB power adapter specifications.

Charge Lights

- The LED will flash white to indicate that the unit is approaching a low battery level.
- The LED will turn red when the battery is depleted, indicating that the unit needs to be charged.
- The LED will glow when charging.
- The LED stays lit when fully charged until the charging cable is unplugged from the device.

Filling the Reservoir

Lift flip top on reservoir and fill the reservoir with lukewarm water.

Note: Always remove charging cable and replace charge port cover prior to filling.




Inserting and Removing Tips



Insert tip at the top of the unit. To remove the tip from the unit, press the tip eject button and pull tip from the unit.

CAUTION: Do not press tip eject button while unit is running.

Adjusting the Pressure Setting

First press power button. Unit will start on "low." Press power button a second time, this will turn unit to "high." Pressing button third time will turn unit OFF. Increase pressure over time to the setting you prefer – or as instructed by your dental professional.

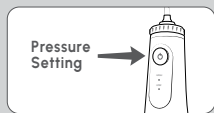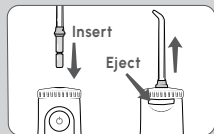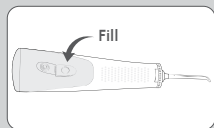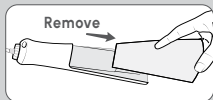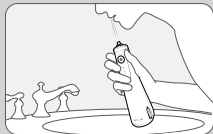Recommended Technique

Lean low over sink and place the tip in mouth. Aim the tip toward teeth, and turn the unit ON. For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth. Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

Tips:

- **LOOK INTO SINK, not in mirror, to reduce splashing.**
- **You can adjust tip direction with a single touch by using the EasySpin™ 360 degree tip rotator.**

When Finished

Turn the unit OFF. Empty any liquid left in the reservoir.

If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove.

Travel mode:

- Locate the power button.
- Press and hold the button down firmly for a continuous 5 seconds until the unit flashes. The device is now in travel mode.
- While in travel mode, the device will not power on.
- To turn off travel mode, press and hold the power button again for a continuous 5 seconds until it flashes.
- After 5 seconds, the device should exit travel mode.
- Verify that the device has returned to normal functionality by checking if any previously disabled features are now accessible.

NOTE: Travel mode will be disabled when unit is plugged in to charge.

NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

SPECIALTY TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets. To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and continue tracing along the gumline.

Orthodontic Tip


To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

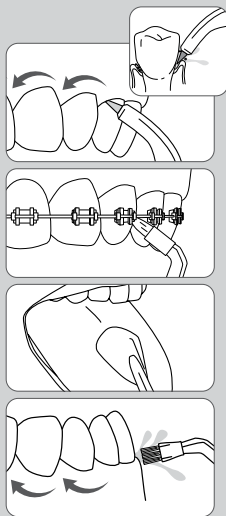Using Mouthwash and Other Solutions

Your Waterpik™ water flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions remove the reservoir (see "When Finished" section) and add from fill door on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

The water reservoir is removable for cleaning and is top rack dishwasher safe. We recommend air drying instead of a heated dry cycle. It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance.

Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON (I) and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water. This process should be done every 1 to 3 months to ensure optimal performance.

Battery Care

If the water flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. Do not expose the battery or appliance to excessive temperatures.

Battery Disposal/Recycling

This product contains a **rechargeable lithium ion battery**. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations.

Please contact your local waste administrative department for any questions.

Note: The rechargeable battery is not replaceable. Remove only at end of product life. This process is not reversible and voids the product warranty.

Battery Removal

To remove the rechargeable battery for recycling, please discharge the battery fully before removing.






Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work.

1. Remove reservoir.
2. Using a Philips screwdriver, remove the screws from the locations shown.
3. Pull the lower cover off.
4. Pull out battery (a) and disconnect from circuit board (b).
5. Recycle/dispose of the battery at your nearest recycling center or take to a hazardous waste collection facility to comply with local and state regulations.

WARNING:

DO NOT cut both battery wires at the same time.

DO NOT disassemble, puncture, hit, squeeze or throw the battery into fire.

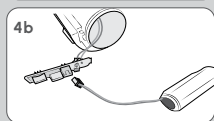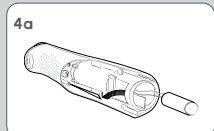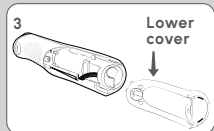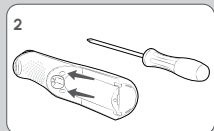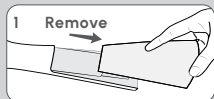Service Maintenance

Waterpik™ water flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs go to www.waterpik.com (US) or www.waterpik.ca (Canada). You may also call our toll-free customer service number at 1-800-525-2774 (US) or 1-888-226-3042 (Canada). Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Still have questions? Go to www.waterpik.com (US) or www.waterpik.ca (Canada).

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

Limitation of Liability

This limited warranty is your sole remedy against Water Pik, Inc. and it shall not be liable for any damages, whether direct, indirect, incidental, special, consequential, exemplary or otherwise, including but not limited to property damage, arising out of any theory of recovery, including statutory, contract or tort.

Arbitration**BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL**

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter. Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim.

Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Colorado and any such arbitration shall be subject to the Federal Arbitration Act and Colorado law.

WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY

Further, by purchasing, installing, or using this product you agree that arbitration must be on an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of Water Pik, Inc. products or the general public. Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

FCC Statement (USA)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by Water Pik, Inc. could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

IC Statement (Canada)

This device complies with Industry Canada's license-exempt RSSs. Operation is subject to the following two conditions: (1) This device may not cause interference; and (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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**Water Pik, Inc., a subsidiary
of Church & Dwight Co., Inc.**